

*"Cutting-edge ECG
technology, when every
heartbeat counts"*

Quark C12-T12X
Stress Testing ECGs



Diagnostic Quality 12-lead Stress Testing ECGs
available in telemetry or cable configuration



COSMED
The Metabolic Company

“*Electrocardiographic (ECG) recording system for continuous monitoring of heart rhythm and evaluation of ischemic ECG changes during exercise and recovery⁽¹⁾*”

- ▶ **Resting and Exercise ECG interpretation**
- ▶ **High resolution ECG processing**
- ▶ **Full Disclosure and Scroll back during the test**
- ▶ **ST segment, ST slopes, Trends**
- ▶ **Retrospective ECG and arrhythmia analysis even during test**
- ▶ **Standard and user defined Exercise Protocols**
- ▶ **Integration with COSMED Gas exchange equipment**
- ▶ **Patient Database & Network Compatibility**



COSMED Stress ECGs are the ideal diagnostic quality PC-based ECG devices for patient care from prevention to rehabilitation. Ultimate quality ECG signal allows an efficient diagnosis of coronary artery disease, the evaluation of patients after artery bypass surgery or angioplasty, and allows to document response to therapies and prevent from myocardial infarction.

Combined with any COSMED gas exchange analysis equipment, COSMED Stress ECGs become the complementary tools that help clinician to quantify patient's functional capacity and evaluate cardiac response in athletes.

COSMED PC-based stress testing systems offer the ultimate in stress testing flexibility, meeting the most advanced requirements for clinical exercise testing.

Two versions are currently available:

- **Quark C12x** - 12-lead Patient Cable Stress Test ECG, with USB communication interface.
- **Quark T12x** - 12-lead Wireless Stress Test ECG (Telemetry). State of the art Digital Radio Frequency technology, immune to environmental interference. Optimized for working in noisy EM Environments (Bluetooth, Wi-Fi Networks, Cellphones etc.). It includes an automatic Channel Shift Feature for selecting the optimal Transmission Channel.

Design

- **12 leads ECG** 12-lead simultaneous recording, display and printout. Detachable 10 lead wires conform to AAMI (AHA) or IEC color codes. Every single lead can be replaced independently
- **High resolution of signal processing** (up to 500 samples/sec) guarantees exceptionally clear tracing for visual view and on screen measurements
- **On-line filtering** Advanced digital filters and automated baseline correction for real-time and stored ECGs to create exemplary ECG traces
- **Integration with COSMED gas exchange equipment** stress ECG and ergospirometry parameters can be thus recorded synchronously, identifying cardiac, pulmonary and vascular function in just a single measurement



(1) AHA Guidelines for Clinical Exercise Testing Laboratories. Circulation. 1995; 91: 912-921




Software Features

Powerful software allows full control over ECG acquisition, real time display, report configuration and printing.

High resolution ECG processing produces an exceptionally clear on-screen display and allows detailed, reliable analysis of ST segments and subtle arrhythmia changes.

Main software features are:

- Full disclosure of entire test enabling post processing on saved data
- Scroll-back during test to see episodes that might have been overseen
- Automatic Arrhythmia detection, print and capture with the possibility to disable unwanted detection types
- Pacemaker detection with mark spikes on screen and in report
- Real time and retrospective J Point and Isoelectric identification
- Predicted oxygen uptake and METS estimation
- Real time super imposition QRST complex
- Remote viewing anywhere on your local hospital network room
- User defined, automatic print during test (standard paper or thermal paper print options)
- Complete reporting functionality allowing the configuration of textual summary and ECG data
- Quality controls including lead quality check and audible signals for lead disconnection and max HR
- User initiated and automatic capture of events
- Pre-programmed user definable testing protocols, including pharmacologic test profiles and many different treadmills and bicycle ergometers
- Export formats: JPEG, XML, GDT, PDF
- Integration with OMNIA Suite including database
- Multiple SpO₂ sensors integration

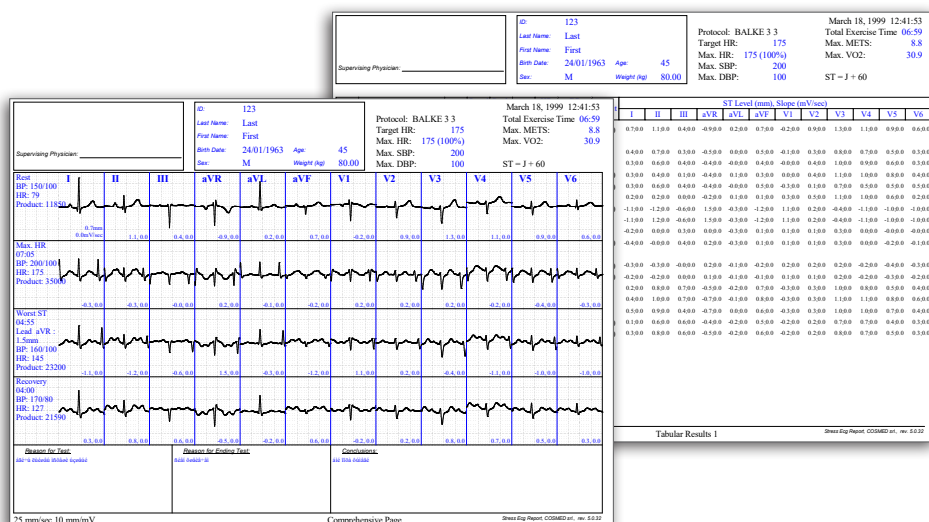

Multiple printout report formats (single lead, all leads, full disclosure, sum up, trends, QRS, etc.)



Nonin® Wrist Ox2 3150 pulse oximeter for SpO₂ integration



Suntech Tango® blood pressure monitor allows to integrate several parameters, including SpO₂

ITALY - Headquarters

COSMED Srl
Rome
+39 06 931-5492
info@cosmed.com

ITALY

COSMED Srl
Milan
+39 02 99765-920
milano@cosmed.com

GERMANY

COSMED Deutschland GmbH
Schweinfurt
+49 (0)9735 81390 00
DE@cosmed.com

FRANCE

COSMED France SASU
Brignais
+33 (0)4 478628053
FR@cosmed.com

THE NETHERLANDS

COSMED Benelux BV
Nieuwegein
+31 (0) 88 10 50 500
BNL@cosmed.com

DENMARK

COSMED Nordic ApS
Odense
+45 6595 9100
DK@cosmed.com

SWITZERLAND

COSMED Switzerland GmbH
Fehraltorf
+41 (0)43 50 869 83
CH@cosmed.com

USA

COSMED USA, Inc.
Concord, Chicago
+1 800 4263763 Toll Free
USA@cosmed.com

AUSTRALIA


COSMED Asia-Pacific Pty Ltd
Artarmon
+61 449 971 170
ANZ@cosmed.com

HONG KONG

COSMED HK Ltd
Kowloon
+852 3708 3126
HK@cosmed.com

CHINA

COSMED Guangzhou Medical
Technology Co. Ltd.
Guangzhou
CHINA@cosmed.com

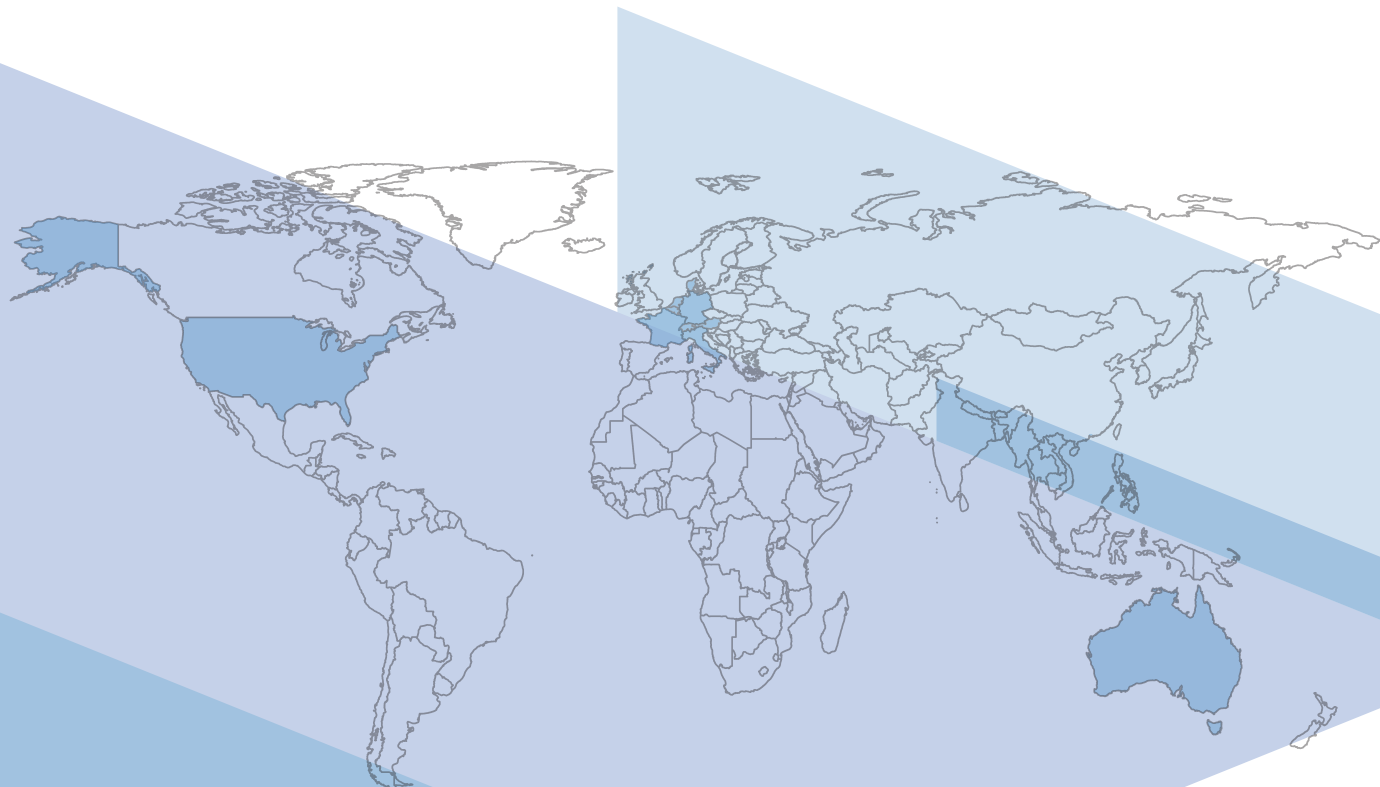


© COSMED

E & OE. Subject to alterations without prior notice.

Products may not be available in your region as depending on countries and certifications.

All trademarks, registered trademarks and logos are the property of their respective owners.

Distributed by



COSMED Srl

Via dei Piani di Monte Savello 37
Albano Laziale - Rome 00041
Italy
+39 (06) 931-5492 Phone
+39 (06) 931-4580 Fax

cosmed.com

To know more:



COSMED ECGs are medical devices.

© 2025/01-a | C02864-02-93